

(12) INTERNATIONAL APPLICATION PUBLISHED UNDER THE PATENT COOPERATION TREATY (PCT)

(19) World Intellectual Property
Organization
International Bureau




(43) International Publication Date
25 March 2004 (25.03.2004)

PCT

(10) International Publication Number
WO 2004/024030 A1

- (51) International Patent Classification⁷: A61F 2/00, A61B 17/00
- (74) Agent: SEVINC, Cenik; Grup Ofis Trademarks & Patents Agency Ltd, Ziya Gokalp Cad. 74/5, 06600 Kolej, ANKARA (TR).
- (21) International Application Number: PCT/TR2003/000014
- (81) Designated State (national): US.
- (22) International Filing Date: 27 February 2003 (27.02.2003)
- (84) Designated States (regional): European patent (AT, BE, BG, CH, CY, CZ, DE, DK, EE, ES, FI, FR, GB, GR, HU, IE, IT, LU, MC, NL, PT, SE, SI, SK, TR).
- (25) Filing Language: English
- (26) Publication Language: English
- (30) Priority Data: 2002/02198 13 September 2002 (13.09.2002) TR
- Published: — with international search report
- (71) Applicant and
(72) Inventor: MALAZGIRT, Zafer [TR/TR]; O.M.U. Department of Surgery, Kurupelit, 55139 SAMSUN (TR).
- For two-letter codes and other abbreviations, refer to the "Guidance Notes on Codes and Abbreviations" appearing at the beginning of each regular issue of the PCT Gazette.

(54) Title: PLUG MADE OF MESH MATERIAL FOR CLOSING LARGE TROCAR WOUNDS

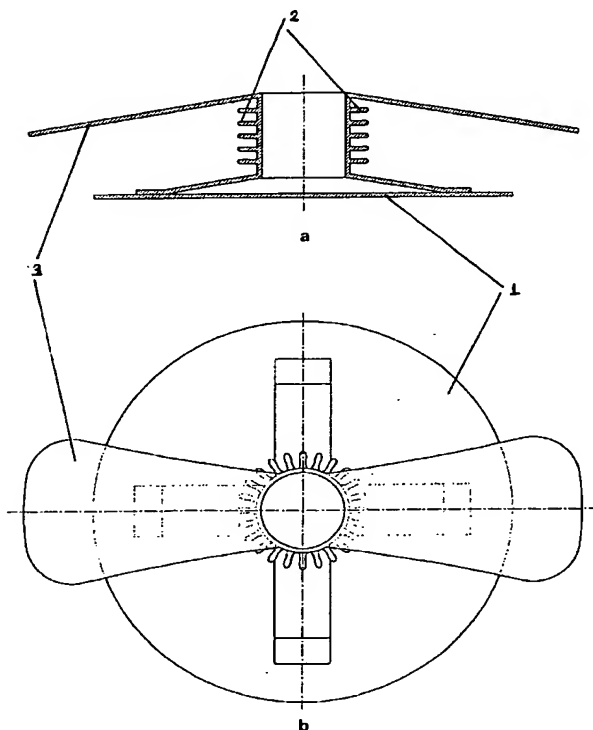

(57) Abstract: The invention is about a material and its relevant technique that is used in the repair of large trocar wounds. The method as a whole is simple to use and minimizes tissue trauma. With this invention, the rest of the suturing techniques are set aside, and by the use of a specially designed plug, the trocar wound repairing process is reduced to a simple and quick procedure. The material is basically made of polypropylene, of which a specific layer covers the intraabdominal face. The material functionally consists of three pieces, and can be called as plug and mesh. The plug and mesh can be used in the repair of almost all large trocar fascial defects that occur at laparoscopic surgery. Due to its unique design, the plug and mesh can fix most of these large trocar wounds safely and securely.

WO 2004/024030 A1

PLUG MADE OF MESH MATERIAL FOR CLOSING LARGE TROCAR WOUNDS

5 TECHNICAL AREA:

The invention is about a plug and mesh material that is to be used in the closure of large trocar wounds without any need for suturing after laparoscopic operations.

10 BACKGROUND OF THE INVENTION:

Trocar site herniation is a recognized complication of laparoscopic surgery. Omental and-sometimes-intestinal herniation with incarceration and obstruction has been documented in recent surgical literature, occurring particularly at 10 mm or larger
15 trocar sites that were not sutured at operation. The necessity to perform fascial closure of any trocar insertion site greater than 5 mm has now been established and is routinely practiced worldwide.

20 However, the conventional closure of such a trocar site fascial defect is often technically difficult, frustrating, indefinitely successful, and even sometimes dangerous due to the limited size of skin incision, the depth of the subcutaneous fatty layer, and necessity of blind manipulation. Moreover, the suturing that
25 involves placement of deep blind sutures after the abdomen has been decompressed is a dangerous manipulation that surgeons tend to avoid.

A number of techniques and instruments have been suggested in the
30 recent 8 to 10 years to facilitate a safe and secure closure of the fascial defect through the tiny skin opening. All of these repairs include passing in any way a suture from one side of the trocar wound to the other, and its ligation. For this purpose

either a heavy needle or a variety of straight needles through which sutures are passed have been used.

Stringer (1995) utilized a Grice needle while Garzotto et al
5 (1995) used a spring-loaded needle. Weiland et al (1995) proposed the use of a Veress needle to pass the sutures through. Robertson (1996) reported a similar method. However, in his method he used two different sutures, of which the ends were taken out through the trocar, and were tied outside, and then pulled outside by
10 pulling one of the sutures. Hellinger et al (1996) have recommended a self-contained disposable fascial closure device. Schwegler et al (1996) described a hook-needle, Thevissen and Pier (1996) reported trocar site closure utilizing the Carter-Thomason technique, and Ammori (1996) recommended the use of a 110-mm heavy
15 spatulated suture needle. In Ammori's method, the ends of a suture, which penetrated the abdominal from skin to skin were brought out of the subcutaneous fat tissue by the help of a nerve hook. Petrakis (1999) recommended a meticulous purse string technique for the closure of trocar defects. Chapman (1999)
20 described the Gore Suture Passer, which are both a needle and a suture grabber.

The basis of the above-mentioned reports can be summarized as that the trocar wound is somehow repaired by any of the suture
25 techniques. Although described in most of the reports as easy and quick methods, these suturing techniques require positioning of the camera and graspers, visualization of the needles during their entrance into the peritoneal cavity, feeding of the graspers or suture passers with the suture loop, all of which are repeated
30 once to thrice for every trocar defect. Any of these suturing techniques are not only time and effort consuming, but also require sophisticated laparoscopic talent and coordination. As more defects at various sites in the abdominal wall are to be

closed after advanced laparoscopic operations, the laparoscopic procedures that support the suturing techniques become more complicated and complex. The above-mentioned suturing techniques would therefore be not that easy and quick.

5

Moreover, a series of manipulation is needed to complete a single suturing. The conventional suturing technique involves much traumatic manipulation including pushing, pulling and retraction of the wound, and insertion and extraction of needles. Most of the time the needle is passed twice, and sometimes more (as depicted in Petrakis' technique). As manipulation in the wound increases, the inflammation and risk of ensuing infection rise considerably. The edema and the collection of seroma and hematoma at the wound further cause dehiscence and hernia formation on a long-term basis.

15

Excessive traumatic manipulation and suturing with heavy sutures oppose the "minimal damage" basis of laparoscopic surgery. The patients are subject to pain and complications at their trocar sites in the postoperative period. The problems associated with the repair of trocar wound would be annoying to the patient as he (or she) is discharged on the first or second postoperative day. The problems of the wound would cause the patient to refer back to the institution.

25

Any of these suturing techniques are to be done under direct vision. It is however impossible to repair the last trocar wound under direct vision. Unless a 0.5 cm scope is used, the last large trocar site can only be closed with conventional blind sutures. At a regular laparoscopic cholecystectomy, the surgeon can only repair the first of two large trocar defects under direct vision. He must close the last one blindly.

30

No matter which suturing technique or needle is used, it is not possible to eliminate the trocar site hernias completely. The current incidence is reportedly around 0.77-3%. As complex laparoscopic surgery becomes more common, the incidence of this complication increases. The reported rates of hernia show that there is not yet any superior method in the safe closure of the trocar fascial defect.

In 2002 Rosin et al advocated closure of the trocar wound by fibrin glue. Their method differs from the previously reported techniques in the following respects. It consists of no sutures, and that it coincides with the "minimally invasive" nature of laparoscopy. The glue can obliterate the subcutaneous tissue and close the skin. However, it seems unrealistic to think that fibrin glue can securely fix the fascial opening.

DETAILED EXPLANATION OF THE INVENTION:

The invention is about a material and its relevant technique that is used in the repair of large trocar wounds. The method as a whole is simple to use and minimizes tissue trauma. With this invention, the rest of the suturing techniques are set aside, and by the use of a specially designed plug, the trocar wound repair is reduced to a simple and quick procedure.

The material is basically made of polypropylene, of which a specific layer covers the intraabdominal face. The material functionally consists of three pieces, and can be called as plug and mesh. This 3-D plug and mesh material agrees well with the rules of hernia surgery, is appropriate to the "minimally invasive surgery" basis of laparoscopy, and is friendly with the physics laws effective in the anterior abdominal wall.

The plug and mesh can be used for the repair of almost all large

trocar fascial defects that occur during laparoscopic surgery. Due to its unique design, the plug and mesh can fix most of these large trocar wounds safely and strongly.

5 Since its use does not necessitate "direct vision", it does not differ whether it is to be used at the first or last trocar wound. Direct vision either from an intraabdominal aspect by laparoscope or from outside by the naked eye is not necessary since the trocar site during insertion of the trocar is evaluated, and relevant
10 preparation of the site is done. To summarize, there is no need to switch to blind suturing during repair of the last trocar wound. The plug and mesh is inserted into place, and its proper positioning and fixation of the ears to the outer fascia are easily accomplished by the guidance of a fingertip.

15

It is easy to use, and its insertion and positioning is very simple and quick. It does not necessitate exhaustive retraction and manipulation. The insertion and placement of the plug in the wound approximately takes a minute or so. If repeated twice in a
20 standard laparoscopic cholecystectomy, or thrice or more in an advanced procedure; it lowers the operating time considerably. The fact that the surgeon repairs the defects easily and safely contributes to a pain-free healing with fewer complications.

25 The plug itself and its insertion and positioning in the wound are natural safety measures against the so-called "chimney effect". The chimney effect is the unwanted escape of the bowel or omentum through the trocar wound during or after the deflation of the pneumoperitoneum. During blind suturing, this effect may cause
30 dangerous complications. The insertion of plug and mesh sweeps away and draws back any herniated intraabdominal structure, and prevents any further herniation thereafter.

It conceals the intraperitoneal orifice (inner ring) and fits perfectly into the trocar's tunnel in the abdominal wall. The specially designed neck of the plug and mesh initially gives a better hold in the wound. The special shaggy design activates in-
5 growth of the granulation tissue, which further stabilizes the plug in place and thus facilitates its better and faster incorporation in the trocar defect.

It is a pressure-friendly material. It functions in company with
10 the physics laws and intraabdominal pressure. In suturing techniques, there is a force-resistance counter-relation between the suture material and tissue tension and pressure. There is no such problem with plug and mesh. The intraabdominal pressure holds the mesh in place, and secures the integrity of the repair. Tissue
15 tension has no detrimental effect on a plug and mesh. On account of Pascal Law, the appropriate surface area of the mesh prevents its prolapse into the defect. Contrary to suturing, the plug and mesh has no physical effects like squeezing, pulling or pushing. This feature contributes to earlier return of the patient to his
20 (or her) normal daily activity, and enables him to even resume sporting activities. This is a fact that conforms well to the basic initiative of the no-tension hernia repair techniques. I specifically recommend its use in patients with increased intraabdominal pressure. Its harmony with Pascal Law makes it a
25 superior choice in patients with debilitating diseases. Its use is especially indicated after laparoscopic surgery in those patients with high intraabdominal pressure (due to ascites, CAPD application in chronic renal failure, and obesity). Another alternative in prevention of hernia formation in those patients
30 with high intraabdominal pressure would be the use of smaller size trocars.

The plug and mesh causes less wound problems in the early

postoperative period. Less tissue trauma decreases the incidence of seroma, hematoma, inflammation, infection and dehiscence, per se. We believe that the figures of hernia incidence would decline by the use of the material, as well.

5

The standard skin stapler or a specially designed stapler can easily fit into the wound through the skin opening. Guided by the fingertip, the ears of the plug and mesh are stapled onto the outer fascial layer one at a time. Since the plug and mesh has two ears, two staples are needed for fixation of each plug, totalling four. Stapling facilitates the repairing by a plug and mesh considerably.

The plug and mesh is made of polypropylene material, a material that causes intraabdominal adhesions if laid in close proximity to intestines or other intraperitoneal structures. The intraabdominal face of the mesh plate should therefore be precovered by an antiadhesive layer. Today there are three composite mesh materials produced by three different companies who have this antiadhesive layer. Namely, these are Parietex Composite (Sofradim International), Sepramesh (Genzyme) and Composix (Bard). All three materials are made of either polypropylene. However, their antiadhesive layers differ. The antiadhesive layer in the first two is absorbable collagen barrier, and ePTFE in the latter.

25

The intraabdominal mesh plate must always lie flat and fully open. It must not be wrinkled or folded up on itself, and it must never come in contact with the intestine in a perpendicular fashion. Such a contact may cause intestinal fistula, and lower its impact on hernia prevention. The special design of the plug and mesh pulls the edges towards the abdominal wall, and prevents the edges from coming in contact with the intestine.

30

Creation of a socket is to be prepared preferably right after the insertion of the trocar if a smooth and clean intraperitoneal surface surrounding the inner trocar hole is not present. At a standard laparoscopic cholecystectomy, such a socket can be formed after insertion of the median epigastric trocar. This can easily be done by a short longitudinal incision along the root of the falciform ligament. The Pneumoperitoneum further blows it open during the procedure, and the socket is ideally formed as the time for insertion of a plug and mesh comes.

The plug and mesh is comprised of three pieces:

1. Intraabdominal mesh plate
2. The shaggy neckpiece
3. The double ear piece.

Intraabdominal mesh plate:

It is made of polypropylene, and its obverse face is covered by either hydrophilic absorbable collagen barrier antiadhesive or ePTFE. It is round in shape.

It has a diameter of 4.5 cm. At calculations we assumed the inner trocar hole to be 1.0 cm² in size. This size will prevent its prolapse into the trocar wound even if the pressure at the trocar site increases by 20 times while it is normal at other points on the mesh. The intraabdominal mesh plate is inserted into the abdominal cavity through the trocar channel, and is laid flat open between the abdominal wall and the intestines centering the inner trocar hole.

Pascal Law says that if we apply a pressure to some part of the surface of a confined fluid by means of a piston, then this pressure will be transmitted without change to all parts of the fluid. Due to the biological properties and consistency of the organs, we can assume that the intraabdominal medium behaves more

or less similar to liquids. Increased intraabdominal pressure will be distributed almost unchanged to all points on the abdominal wall.

5 According to the Pascal Law, the pressures at two pistons are equal. However, the force is directly correlated with the surface area of the piston. Normally, the pressure on the inner trocar hole is similar with the pressure at different points on the mesh, a fact that prevents the herniation of the mesh into the trocar
10 channel.

Let us consider now an odd condition that the pressure on the inner trocar hole is increased to a peak of 20 times that of other points on the mesh. In order to prevent prolapsing of the mesh
15 into the trocar tunnel, the total force on the mesh must be higher than or equal to the force on the inner trocar hole.

$$\Delta P = \Delta F / \Delta A$$

ΔP : pressure at a certain point

$$\Delta P_{th} = \Delta P_{ps}$$

ΔF : force at a certain point

$$20 \quad \Delta F_{th} / \Delta A_{th} = \Delta F_{ps} / \Delta A_{ps}$$

ΔA : area influenced by the force

$$\Delta F_{ps} = 20 \cdot \Delta F_{th}$$

th: trocar hole

$$\Delta F_{th} / \Delta A_{th} = 20 \Delta F_{th} / \Delta A_{ps}$$

ps: plate surface

$$\Delta A_{ps} = 20 \cdot \Delta A_{th}$$

Since the mesh is small in size, πr^2 will be used in the
25 calculation.

$$\pi r_{ps}^2 = 20 \cdot \pi r_{th}^2$$

$$r_{ps}^2 = 20 \cdot r_{th}^2 \quad (r_{th} = 0.5 \text{ cm})$$

$$r_{ps}^2 = 20 \times 0.25 \text{ cm}^2$$

$$30 \quad r_{ps}^2 = 5.0 \text{ cm}^2$$

$$r_{ps} = 2.24 \text{ cm} \cong 2.25 \text{ cm}$$

$$R_{ps} = 4.5 \text{ cm.}$$

This final figure shows what diameter the plate should be. For bigger trocar holes, the size of the plate must be recalculated accordingly.

5 The shaggy neck-piece:

This piece, which is made of polypropylene connects the plate to the double-ear piece. Its specially designed shaggy structure helps in taking a good hold of the trocar tunnel.

- 10 It is a cylinder with a diameter and height of 1cm. The shaggy appearance comes from the projecting loops of polypropylene very similar to those of a heavy towel. The loops are 2 mm in length. Thus the diameter of the neck reaches to 14 mm in total. Towards bottom the neck cylinder opens up to form four legs, all of which
15 unite with the plate's reverse side. The shaggy neckpiece joins with the two ears above.

The bottom of the cylinder extends to form four legs at an angle of 100°. Every foot is 1.3 cm in length and 0.6 cm in width. All
20 the feet merge with the plate at its periphery with 0.3 cm segments at the tip. When the ears are pulled upwards too much or too strongly, the nails and feet pull up the plate from its periphery. Thus the ends of the plate are collected upwards, inwards and away from the intestine. If these feet are not
25 existing, and the neck pulls the plate directly from its center, then the plate will be pulled upwards bending the ends down towards the intestine. This causes perpendicular contact of the edges with the intestine. This is the most unwanted position of the mesh.

30

Above, the shaggy neck cylinder opens up to both sides to form the ears, which look very much like a plane propeller. The circumference of the neck cylinder is 3.14 cm. The radix of each

ear is 1.0 cm in width. There are two 0.7 cm bare neck edges between the ears.

Double-ear piece:

- 5 The ears extend to the opposite directions. Each ear is 1.0 cm in width and 3.0 cm in length from the center of the neck cylinder. It reaches to a width of 1.7 cm at the end of this distance. It starts from the neck with an angle of 90°. It bends downwards to form an overall angle of 70°. At its radix its shape is concave,
10 but becomes flat laterally. This concavity and the parabolic angling give the material an extra holding strength in the trocar tunnel. The ears function like a lever, and pull the plate up to a proper tightness. The neck and ears contribute to the general behavior of the material as if it is a worker coming out of an
15 underground canal in the street.

How is the material (invention) used?

- At the end of any laparoscopic operation, the surgeon decides on which trocar wounds are to be repaired by the plug and mesh. The
20 surgical nurse should prepare the equivalent number of plug and mesh material and the "multipurpose stapler" (if available) at her table.

- Using the laparoscope, the surgeon rechecks the inner hole whether
25 it has a clean flat area around with a radius of 2.5 cm. Any hole, around which a sufficient clean flat area is impossible to develop, must be excluded from the plug and mesh repair.

- It is easier to place a plug and mesh material in the presence of
30 pneumoperitoneum, however this is not essential. It can be used safely after complete deflation, as well. If the surgeon wants to keep the pneumoperitoneum before and after insertion of the first plug, one of his (or her) assistants must employ a finger (Figures

2b, 2c). If not done so, some gas will eventually escape through the mesh.

Insertion of a Plug and Mesh does not require any form of direct vision. However, inner vision of the insertion and placement of the mesh plate can be visualized by the laparoscope if it is in the abdomen and not used for any other purpose. In the first few cases, the surgeon may not exactly know with how much force he should pull the ears up. It is helpful to the surgeon in this initial phase for developing a sensation of an ideal pullback tension. After the first few cases, insertion and positioning of the material can easily be accomplished blindly. The last trocar wound to be repaired must preferably be the umbilical one.

The subcutaneous tissue above the outer fascial layer adjacent to the trocar tunnel is slightly dissected by a middle-sized clamp to form two small gaps (Figures 3a, 3b). During this step the size of the gaps must be memorized for the next step, at which point the surgeon must decide whether he should trim the ears. The direction of the gaps must be parallel to the fascial tension lines. In the anterior abdominal wall, the gaps are placed horizontally.

A heavy curved clamp holds the plug and mesh in such a position that the ears are on top of each other (Fig 4a). The tip of the clamp is exactly at the neck-ear junction. The nurse soaks the material into saline solution. As the assistant pulls away his finger, the surgeon inserts the clamped material into the wound, and gently pushes it with a steady jerk until the bulk gets into the abdominal cavity (Fig 4b). The resistance suddenly disappears as its entire length enters into the abdominal cavity. Pushing it further to various directions or to-and-fro movements in the abdomen is not recommended. At this point the surgeon pulls back the clamp until he sees in the screen or feels at his hand that

the mesh plate is touching the parietal peritoneum slightly. After checking that the tips of the ears are out of the fascial layer, the clamp is taken out. The ears are brought all the way out of the skin (Figure 6a).

5

A clamp checks the gaps that were prepared beforehand. The ears are trimmed accordingly. The trimmed tip of an ear is held by a middle-sized clamp, and pushed inside and to the gap settling onlay on the fascia (Figure 6b). The surgeon's index fingertip helps accomplishing this step very easily. The same is repeated for every trocar site.

If a "multi-purpose stapler" is available, it is inserted into the gap, and is shot to put down a titanium clip over the mesh to the fascia. Ideally, the clip should be put at outer one-third segment of each ear.

As the ear is put into its gap, the fingertip checks its proper position. With the finger in position, the stapler may be pushed in and fired. If both cannot fit into the hole, then stapler is pushed while the fingertip is pulled back. The stapler is fired halfway to let the clip appear at the nozzle. The stapler then tackles the mesh, and both are stretched laterally, and the stapler is fully fired. The same is repeated on the other side.

25

SHORT DESCRIPTIONS OF THE FIGURES:

Figure 1a : Side view of the plug and mesh.

Figure 1b : View of the plug and mesh from above.

Figures 2a, 2b, 2c : Extraction of the trocar at the end of the primary procedure, and insertion of the finger into the wound to halt escape of the gas.

Figures 3a, 3b : The making of a gap in the subcutaneous tissue, in which thereafter the double ear piece (3) will be settled in.

Figures 4a, 4b : The insertion of the mesh and plug in the trocar hole. Note that the clamp holds the material at its double ear segment (3).

5 Figures 5a, 5b : The positioning of the mesh plate (1) over the inner trocar hole by pulling the double ear piece upwards.

Figures 6a, 6b, 6c : The placement of the double ear piece (3) into the subcutaneous tissue, and its fixation.

Figure 7 : The final appearance of the plug and mesh in the trocar hole.

10

CLAIMS:

1. It is the plug and mesh that is to be used for the repair of large trocar wounds; it is formed from three pieces that;

5 a. is the intraabdominal mesh plate (1) that it is made of polypropylene, and that its obverse face is covered by either hydrophilic absorbable collagen barrier antiadhesive or ePTFE, that it is round in shape, and that it is laid between the abdominal wall and the intestines centering the inner trocar hole.

10 b. is the shaggy neck piece (2) that it is made of polypropylene, that it connects the intraabdominal mesh plate (1) to the double-ear piece (3), that it has a shaggy or thorny outer surface, that it fits into the trocar hole, that it is cylindrical in shape, and that it gets a better hold
15 due to its shaggy shape.

c. is the double-ear piece (3) that it is made of polypropylene, that its two ears extend in opposite directions from the shaggy neck piece (2), that it has a shape of a propeller, and that it is fixated to the outer layer of the fascia.
20

2. It is the shaggy neck piece (2) that it is connected to Claim 1, that its bottom opens up at a certain angle to form four legs, all of which unite with the intraabdominal mesh plate.

25 3. It is the double ear piece in connection with Claim 1 that it has a concave shape at its radix, but becomes flat laterally.

4. It is the double-ear piece in connection with Claim 1 and 3 that it has a parabolic angling at its vertical section, and that
30 every ear of it has a certain angle with the cylindrical body.

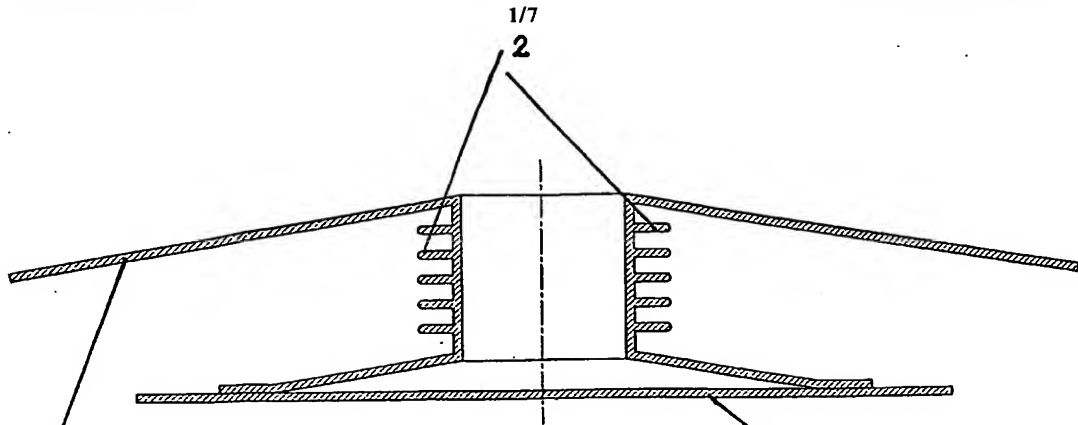


FIGURE 1a

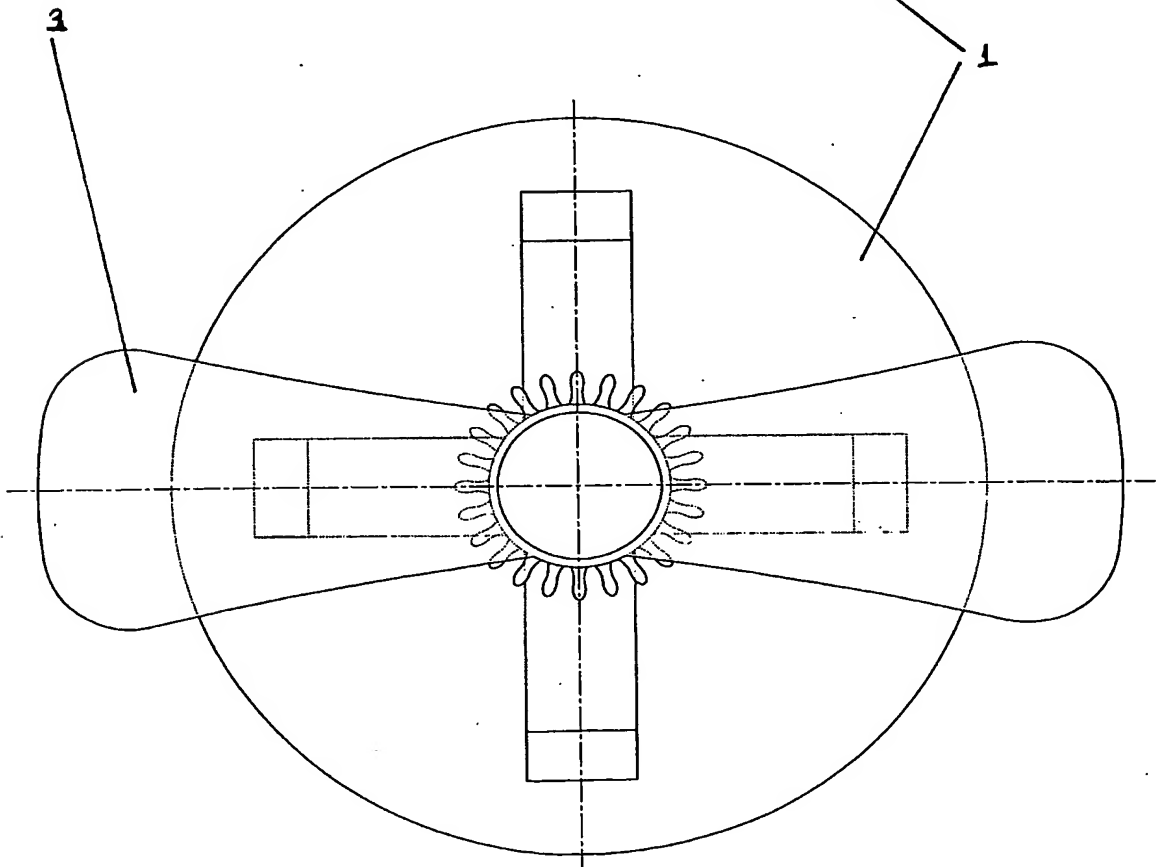


FIGURE 1b

2/7

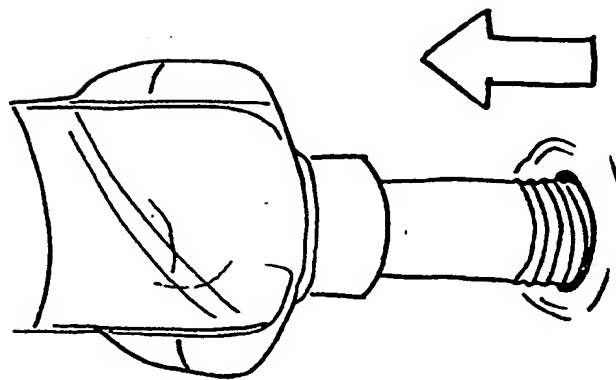


FIGURE 2a

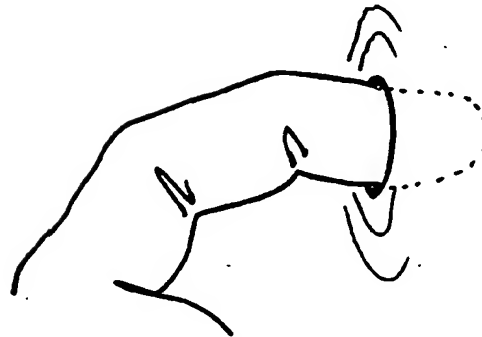


FIGURE 2b

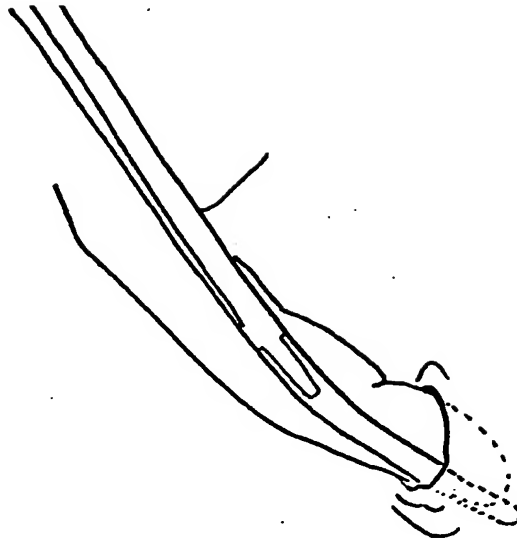


FIGURE 2c

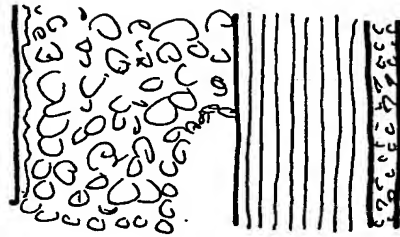


FIGURE 3b

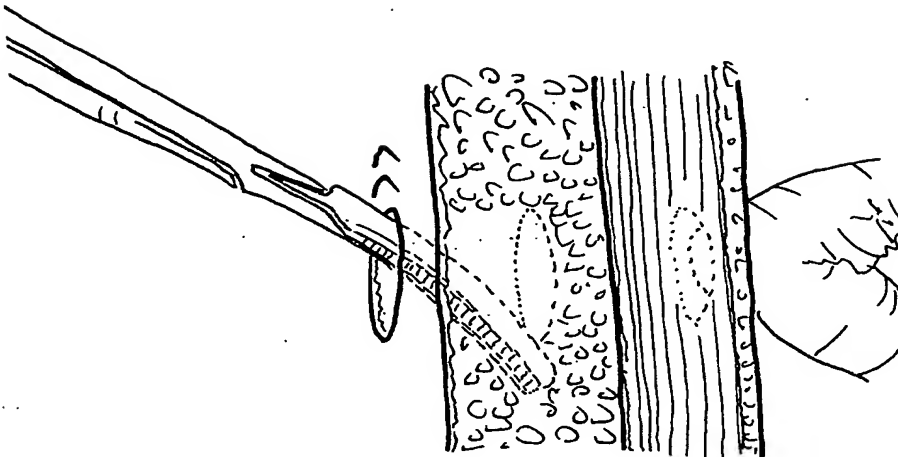


FIGURE 3a

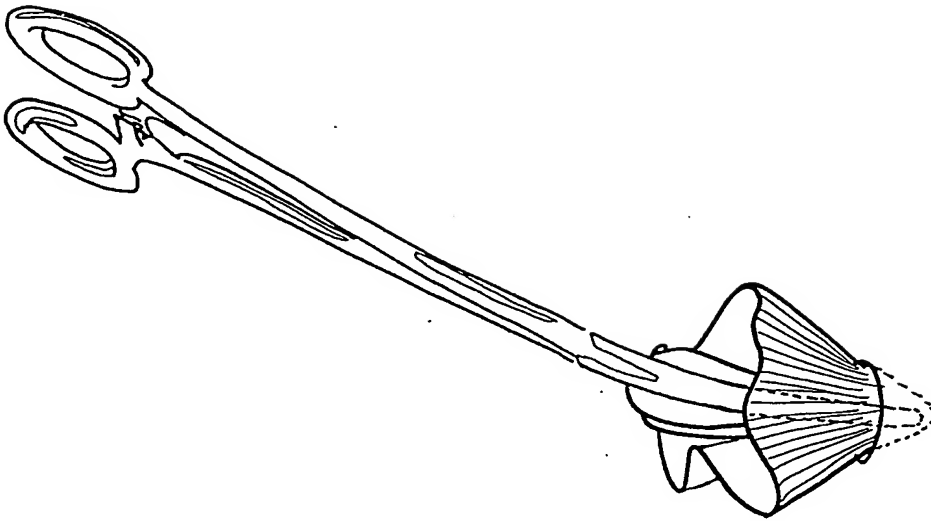


FIGURE 4b

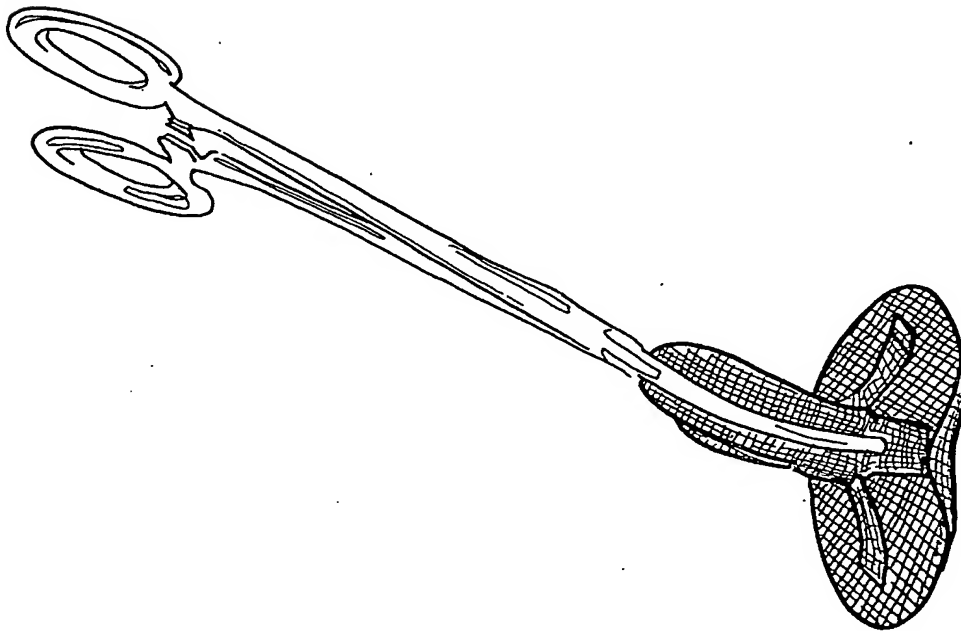


FIGURE 4a

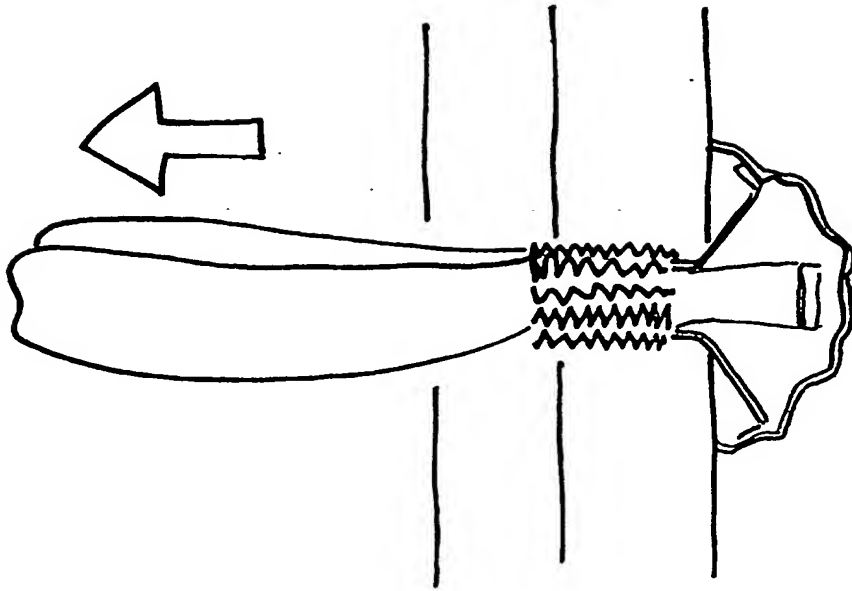


FIGURE 5a

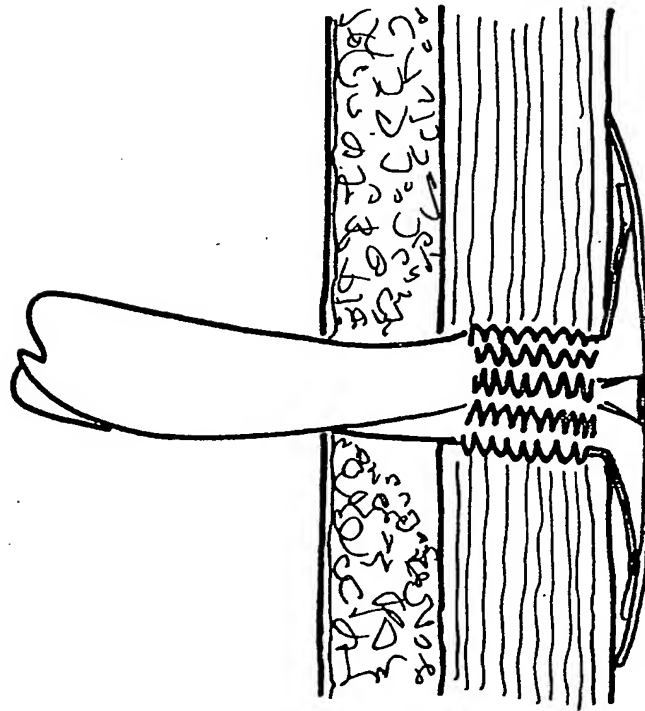


FIGURE 5b

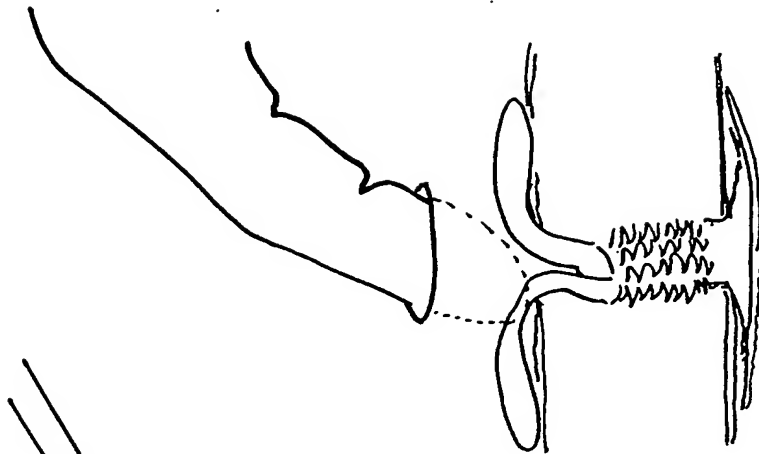


FIGURE 6b

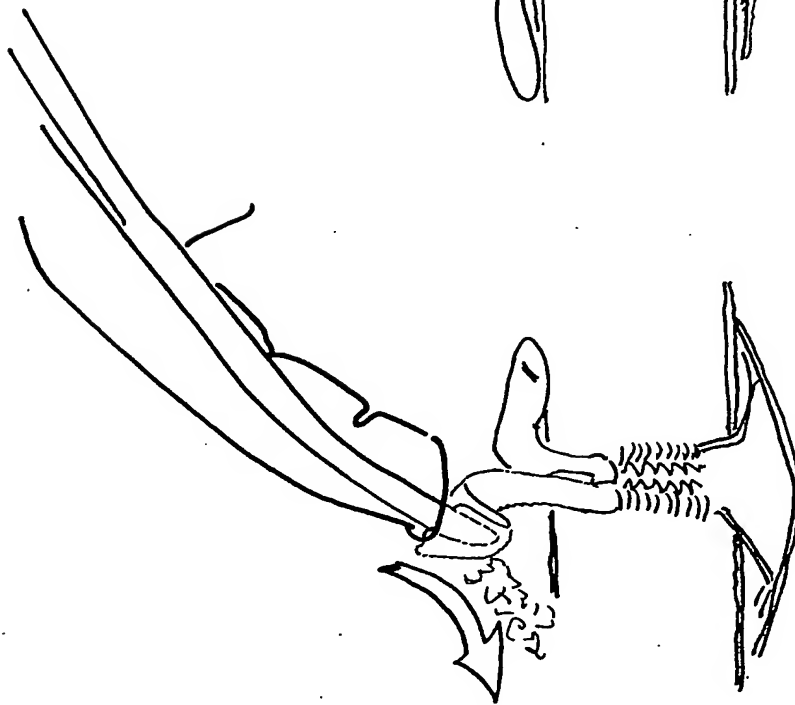


FIGURE 6a

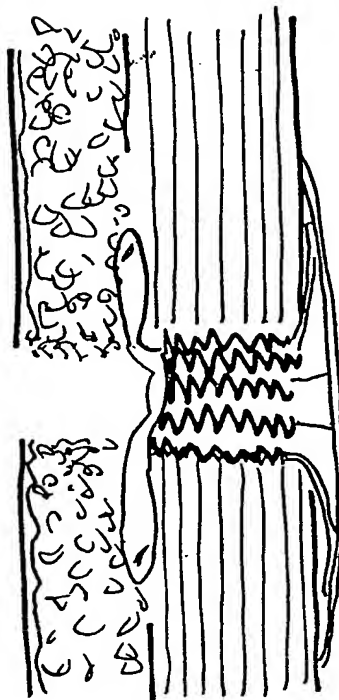


FIGURE 7

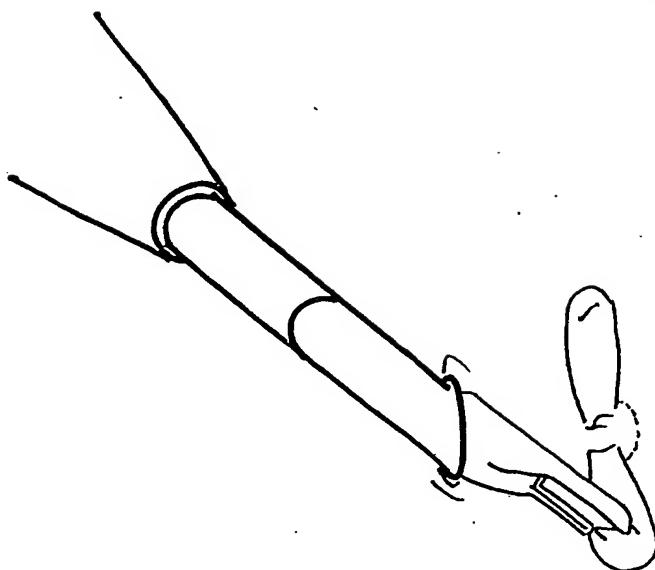


FIGURE 6c

INTERNATIONAL SEARCH REPORT

International Filing No
PCT/TR 03/00014

A. CLASSIFICATION OF SUBJECT MATTER
IPC 7 A61F2/00 A61B17/00

According to International Patent Classification (IPC) or to both national classification and IPC

B. FIELDS SEARCHED

Minimum documentation searched (classification system followed by classification symbols)
IPC 7 A61F A61B

Documentation searched other than minimum documentation to the extent that such documents are included in the fields searched

Electronic data base consulted during the international search (name of data base and, where practical, search terms used)

EPO-Internal

C. DOCUMENTS CONSIDERED TO BE RELEVANT

Category *	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.
X	EP 0 898 944 A (ETHICON INC) 3 March 1999 (1999-03-03) column 3, line 15 - line 37; figures 1,5	1
A	WO 00 74576 A (BURGARD GUNTHER ;NOVOMED GMBH (DE)) 14 December 2000 (2000-12-14) page 8, paragraph 6 -page 9, paragraph 4; figures 2,3	1
A	US 5 304 117 A (WILK PETER J) 19 April 1994 (1994-04-19) column 4, line 42 -column 5, line 27; figures 4B,5	1
A	US 2001/027347 A1 (ROUSSEAU ROBERT A) 4 October 2001 (2001-10-04) paragraph '0049!; figure 9	1

☐ Further documents are listed in the continuation of box C.

☒ Patent family members are listed in annex.

* Special categories of cited documents :

A document defining the general state of the art which is not considered to be of particular relevance

E earlier document but published on or after the international filing date

L document which may throw doubts on priority claim(s) or which is cited to establish the publication date of another citation or other special reason (as specified)

O document referring to an oral disclosure, use, exhibition or other means

P document published prior to the international filing date but later than the priority date claimed

T later document published after the international filing date or priority date and not in conflict with the application but cited to understand the principle or theory underlying the invention

X document of particular relevance; the claimed invention cannot be considered novel or cannot be considered to involve an inventive step when the document is taken alone

Y document of particular relevance; the claimed invention cannot be considered to involve an inventive step when the document is combined with one or more other such documents, such combination being obvious to a person skilled in the art.

& document member of the same patent family

Date of the actual completion of the international search

18 September 2003

Date of mailing of the international search report

06/10/2003

Name and mailing address of the ISA
European Patent Office, P.B. 5818 Patentlaan 2
NL - 2280 HV Rijswijk
Tel (+31-70) 340-2040, Tx. 31 651 epo nl,
Fax (+31-70) 340-3016

Authorized officer

Angeli, M

INTERNATIONAL SEARCH REPORT

Information on patent family members

International / tion No
PCT/TR 03/00014

Patent document cited in search report		Publication date	Patent family member(s)	Publication date
EP 0898944	A	03-03-1999	US 6241768 B1 EP 0898944 A2 JP 11146888 A	05-06-2001 03-03-1999 02-06-1999
WO 0074576	A	14-12-2000	DE 29909888 U1 AU 5812600 A WO 0074576 A1 EP 1182973 A1	19-10-2000 28-12-2000 14-12-2000 06-03-2002
US 5304117	A	19-04-1994	NONE	
US 2001027347	A1	04-10-2001	US 6425924 B1 CA 2342641 A1 EP 1219265 A2 JP 2002209901 A CA 2342640 A1 EP 1145693 A2 JP 2001309933 A US 2002147457 A1	30-07-2002 29-06-2002 03-07-2002 30-07-2002 30-09-2001 17-10-2001 06-11-2001 10-10-2002